


AquaTop XS

Electronic warm water irrigator



The AquaTop XS is an ear irrigator available in two versions. One is suitable for ear cleaning purposes with its water stream at a constant 37 degrees. The other is a three temperature version for ear cleaning and caloric simulation. This has three temperature settings (30°C, 37°C and 44°C). The electronic temperature control is highly accurate. The AquaTop XS incorporates a constant flow design to keep fresh water flowing through the irrigator up to the handpiece to prevent biofilm and microbe growth, making it safe for doctor and patient.

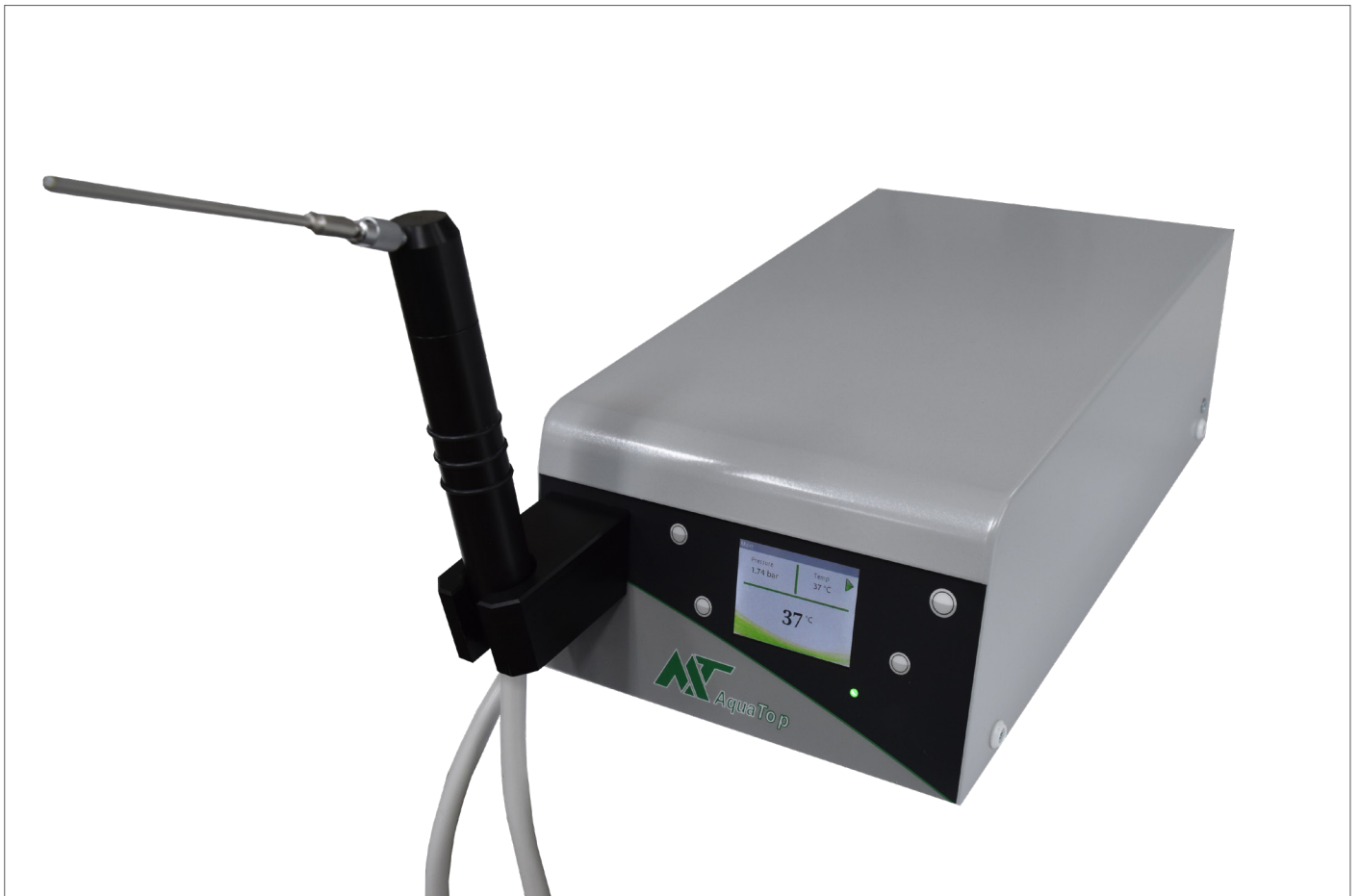

Benefits at a glance

- Constant water flow system to prevent microbe growth and biofilm
- Accurate electronic temperature control
- Pressure regulation from button on handpiece
- Overheating protection
- LCD display
- LED indicates device status
- High quality construction
- Easy integration with the MediTop ENT treatment units
- Compact design



Technical specifications

Voltage	230 Volt AC - 50/60 Hz
Temperature settings	30°C - 37°C - 44°C
Minimum inlet water pressure	2.0 bar
Maximum inlet water pressure	6.0 bar
Maximum power consumption	1100 VA
Water pressure setting	1.8 bar
Weight	3 kg