

PATIENTSTOL 5104 FRÅN TYSKA JÖRG UND SOHN





- ♦ Universal patient chair for ENT and other medical examinations.
- ♦ In modern design.
- ♦ With infinitely variable height-adjustment by means of a double pedal switch.
- ♦ Maintenance-free, electromotive lift, stroke 200 mm with 12-V-motion according IEC 601-1.
- ♦ Horizontally installed seat (slight seat inclination).
- ♦ The backrest adjustable from the vertical (95°) up to and below the horizontal position by gaspring.
- ♦ Infinitely variable height adjustment of the headrest. Can be completely pulled out if necessary.
After loosening the hand screw additionally to be adjusted to approx. 170° incline.
- ♦ The complete legrest folding up to horizontal position. Legrest extensible by synchronously following-on footrest.
- ♦ The two comfortably upholstered armrests are tiltable forwards and backwards.
- ♦ The complete upper part fully rotating and progressively lockable.
- ♦ The one-piece upholstery is fixed with Velcro fastener.
- ♦ Option: zeroizing pedal switch.

Technical datas:

Height adjustment: 510 – 710 mm

Voltage: 230 V, 50 cycles

Current consumption: 2,7 A

Weight: approx. 80 kg


Subject to technical modifications!



Patient chair

5104

Dimensioned drawing:

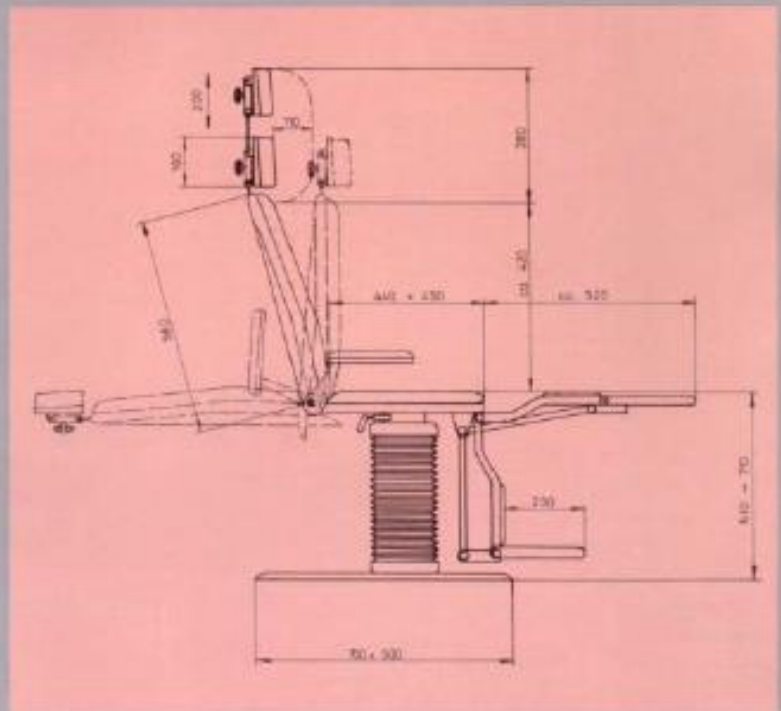

jörg&sohn

Dental
Equipment

Medical Equipment

Podiatry

Lindenstraße 24 · 96482 Ahorn-Schorkendorf · Tel. 095 61/82 99-0 · Fax 095 61/82 99-37

Philip Ströberg
OPIVA Nordic
070-971 81 40
philip.stroberg@opivanordic.se

Mette Ståhlgren
OPIVA Nordic
0725-30 28 82
mette.stahlgren@opivanordic.se

Niclas Holmebäck
OPIVA Nordic
0725-29 28 82
niclas.holmeback@opivanordic.se